


Bicarb System

The AmeriWater Bicarb System utilizes virgin polyethylene tanks with cone bottom, complete drain fittings, and gasketed 12" covered fill hole. Piping and true union valves are PVC schedule 80. Pumps are 115V, single-phase 8.5 amp motors, with all stainless steel design. The control is 115V and has timed fill for the mix tank, timed mixing and filling control with low and high level alarms. The high vortex pumped mixing action provides complete mixing in minutes and keeps the bicarbonate powder in solution.

Lightly pressurized distribution and recirculation provides steady, consistent delivery and facilitates complete and efficient disinfection.

Features:

- Automatic fill and mix
- Lightly pressurized for complete recirculation of the distribution loop
- Complete, quick, and easy disinfection
- Inlet flow restrictor and manual by-pass on pure water filling
- Ozone disinfection connections

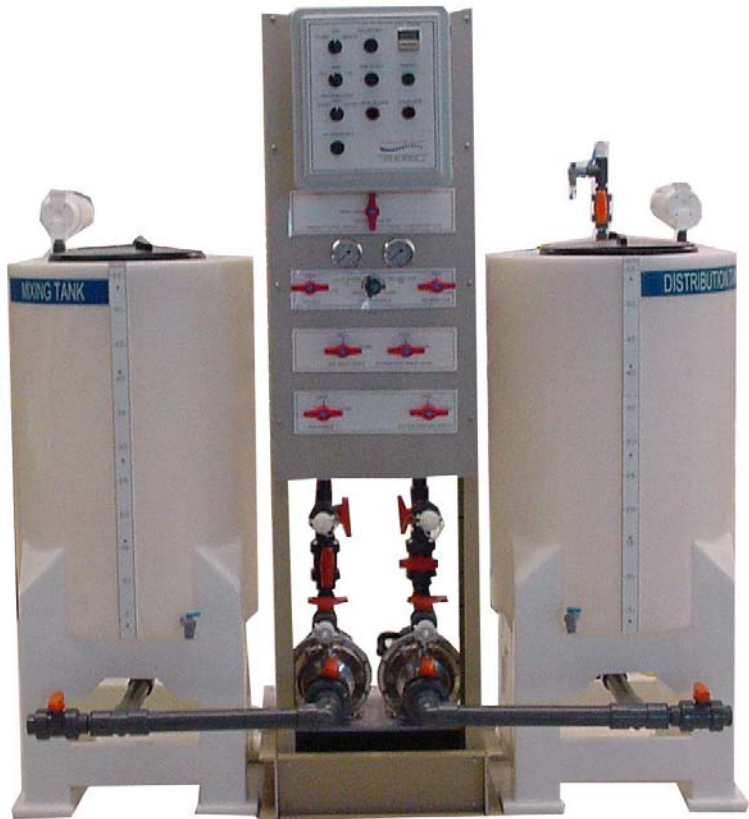

Bicarb Mix and Distribution System Dual Tank System

Part Number	Description
00BC55-55	55 Gallon Mix and Distribution
00BC100-100	100 Gallon Mix and Distribution
00BC100-200	100 Gallon Mix/200 Gallon Distribution

Bicarb Mix Single Tank System

Part Number	Description
00BC100	100 Gallon Single Tank Mixing with Easy Jug-Filling Port or Pumped Distribution